

ELISA Flex: Mouse IFN- α (HRP)

3326-1H-6 | 3326-1H-20

ELISA Flex kit for quantitative determination of native and recombinant mouse IFN- α in solution, e.g. cell supernatant and serum/plasma.

The kit includes		3326-1H-6 for 6 plates	3326-1H-20 for 20 plates
Capture mAbs:	MT24A/44A/104 (0.5 mg/ml)	350 μ l	1200 μ l
Detection mAbs:	MT9L/14A/113, biotinylated (0.5 mg/ml)	150 μ l	500 μ l
Streptavidin-HRP		80 μ l	250 μ l
Recombinant mouse IFN- α ELISA standard		1 vial	1 vial
Standard reconstitution buffer A5		1 ml	1 ml


To ensure total recovery of the stated quantity, vials have been overfilled.

Shipping and storage

Shipped at ambient temperature. All reagents should be stored at 4-8 °C upon receipt, except the standard which should be stored at -20 °C. Antibodies are supplied in sterile-filtered PBS with sodium azide (0.02%). Streptavidin-HRP is supplied in PBS with 0.002% Kathon CG. The expiry date indicates how long unopened products, stored according to instructions, are recommended for use.

Protocol

- Day 1**
1. Add 100 µl/well of capture mAbs MT24A/44A/104 diluted to 3 µg/ml in PBS, pH 7.4. Use high protein binding ELISA plates. Incubate overnight at 4–8 °C.
- Day 2**
2. Empty the plate and add 200 µl/well of PBS with 0.05% Tween 20 and 0.1% BSA (incubation buffer) to block the plate. Incubate for 1 hour at room temperature.
 3. Wash the plate 5 times with PBS containing 0.05% Tween 20 (300 µl/well).
 4. Add 100 µl/well of samples or standards diluted in incubation buffer. Include assay background control, i.e. wells without standard. Incubate for 2 hours at room temperature.
 5. Wash as above.
 6. Add 100 µl/well of detection mAb MT9L/14A/113-biotin diluted to 1 µg/ml in incubation buffer. Incubate for 1 hour at room temperature.
 7. Wash as above.
 8. Add 100 µl/well of Streptavidin-HRP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature. Please note that sodium azide used in buffers will inhibit HRP activity.
 9. Wash as above.
 10. Add 100 µl/well of TMB substrate (product code: 3652-F10) and incubate for 15 minutes.
 11. Add 100 µl/well of 0.2 M H₂SO₄ to stop the reaction.
 12. Measure the optical density in an ELISA reader at 450 nm within 15 min. Preferably use a reader capable of subtracting a reference wavelength of between 570 and 650 nm. Representative standard curve shown below.

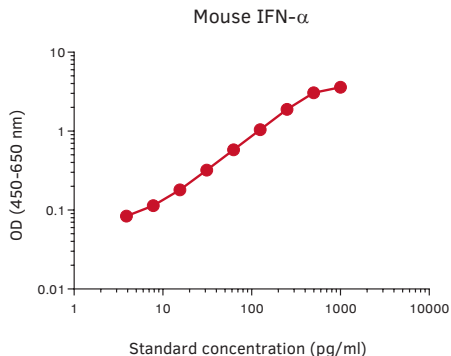

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.



The products are for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office)
Sweden
Tel: +46 8 716 27 00
mabtech@mabtech.com

Mabtech, Inc.
USA
Tel: +1 513 871-4500
mabtech.usa@mabtech.com