

## ELISA Flex: Porcine IL-5 (ALP)

3194-1A-6 | 3194-1A-20

ELISA Flex kit for quantitative determination of native and recombinant porcine IL-5 in solution, e.g. cell supernatant.

The kit includes		3194-1A-6 for 6 plates	3194-1A-20 for 20 plates
Capture mAb:	MTAL25 (0.5 mg/ml)	300 µl	1000 µl
Detection mAb:	MTAL59, biotin (0.5 mg/ml)	150 µl	500 µl
Streptavidin-ALP		80 µl	250 µl
Recombinant porcine IL-5 ELISA standard		1 vial	1 vial
Standard reconstitution buffer A8		1 ml	1 ml

To ensure total recovery of the stated quantity, vials have been overfilled.

### Shipping and storage

Shipped at ambient temperature. All reagents should be stored at 4-8 °C upon receipt, except the standard which should be stored at -20 °C. Antibodies are supplied in sterile-filtered PBS with sodium azide (0.02%). Streptavidin-ALP is supplied in 0.1 M Tris buffer with 0.002% Kathon CG. The expiry date indicates how long unopened products, stored according to instructions, are recommended for use.

# General and Preparations

## Specificity

The kit contains a matched pair of monoclonal antibodies (mAbs) specific for native and recombinant porcine IL-5.

## Standard range

10-2000 pg/ml

## Calibration


No international standard exists for calibration.

## Reconstitution of ELISA standard

Reconstitute the ELISA standard to a stock solution of 0.5 µg/ml by adding 1 ml of the standard reconstitution buffer. Allow the standard to dissolve for 5 minutes and mix thoroughly. The standard should be kept in aliquots at -20 °C. Avoid repeated freeze-thaw cycles.

## Preparation of standard curve

Prepare within 30 minutes of use. Volumes are sufficient for duplicates.

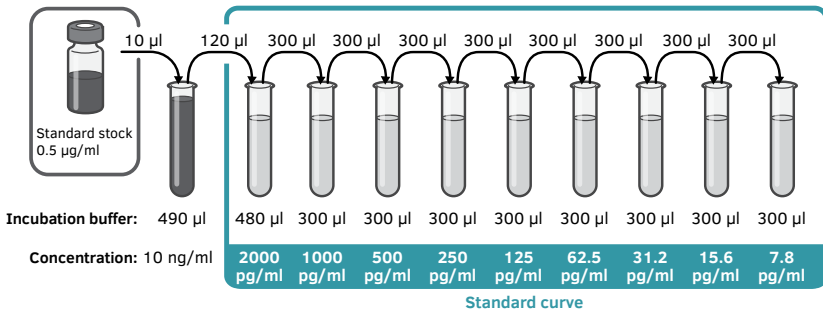

# Protocol

## Day 1

1. Add 100  $\mu$ l/well of capture mAb MTAL25 diluted to 2  $\mu$ g/ml in PBS, pH 7.4. Use high protein binding ELISA plates. Incubate overnight at 4-8 °C.

## Day 2

2. Empty the plate and add 200  $\mu$ l/well of PBS with 0.05% Tween 20 and 0.1% BSA (incubation buffer) to block the plate. Incubate for 1 hour at room temperature.
3. Wash the plate 5 times with PBS containing 0.05% Tween 20 (300  $\mu$ l/well).
4. Add 100  $\mu$ l/well of samples or standards diluted in incubation buffer. Include assay background control, i.e. wells without standard. Incubate for 2 hours at room temperature.
5. Wash as above.
6. Add 100  $\mu$ l/well of detection mAb MTAL59-biotin diluted to 1  $\mu$ g/ml in incubation buffer. Incubate for 1 hour at room temperature.
7. Wash as above.
8. Add 100  $\mu$ l/well of Streptavidin-ALP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature.
9. Wash as above.
10. Add 100  $\mu$ l/well of pNPP substrate (product code: 3652-P10) and incubate at room temperature protected from direct light for approximately 60 minutes.
11. Measure the optical density in an ELISA reader at 405 nm. Preferably use a reader capable of subtracting a reference wavelength of between 570 and 650 nm. Representative standard curve shown below.

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.



**The products are for research use only.**

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

**Mabtech AB (Head Office)**  
Sweden  
Tel: +46 8 716 27 00  
mabtech@mabtech.com

**Mabtech, Inc.**  
USA  
Tel: +1 513 871-4500  
mabtech.usa@mabtech.com