

ELISA Flex: Human IFN- γ (HRP)

3420-1H-6 | 3420-1H-20

ELISA Flex kit for quantitative determination of native and recombinant human IFN- γ in solution, e.g. cell supernatant and serum/plasma.

The kit includes		3420-1H-6 for 6 plates	3420-1H-20 for 20 plates
Capture mAb:	1-D1K (1 mg/ml)	150 μ l	500 μ l
Detection mAb:	7-B6-1, biotin (1 mg/ml)	80 μ l	250 μ l
Streptavidin-HRP		80 μ l	250 μ l
Recombinant human IFN- γ ELISA standard		1 vial	1 vial
Standard reconstitution buffer A8		1 ml	1 ml

To ensure total recovery of the stated quantity, vials have been overfilled.

Shipping and storage

Shipped at ambient temperature. All reagents should be stored at 4-8 °C upon receipt. Antibodies are supplied in sterile-filtered PBS with sodium azide (0.02%). Streptavidin-HRP is supplied in PBS with 0.002% Kathon CG. The expiry date indicates how long unopened products, stored according to instructions, are recommended for use.

General and Preparations

Specificity

The kit contains a matched pair of monoclonal antibodies (mAbs) specific for native and recombinant human IFN- γ .

Standard range

4-400 pg/ml

Calibration

The ELISA standard has been calibrated against an international standard from the National Institute of Allergy and Infectious Diseases (NIAID) Bethesda, US. One ng of supplied standard equals 39 U Gxg01-902-535 NIAID-standard. Please note that the calibration is batch specific.

Analysis of serum and plasma samples


Analysis of serum/plasma requires the use of ELISA diluent (product code: 3652-D2). The ELISA diluent blocks heterophilic antibodies, commonly found in serum/plasma, from cross-linking the assay antibodies, thereby preventing false positive read-outs. The ELISA diluent should be used for dilution of standard, samples, and detection antibody.

Reconstitution of ELISA standard

Reconstitute the ELISA standard to a stock solution of 0.2 $\mu\text{g/ml}$ by adding 1 ml of the standard reconstitution buffer. Allow the standard to dissolve for 5 minutes and mix thoroughly. The standard should be kept in aliquots at $-20\text{ }^{\circ}\text{C}$. Avoid repeated freeze-thaw cycles.

Preparation of standard curve

Prepare within 30 minutes of use. Volumes are sufficient for duplicates.

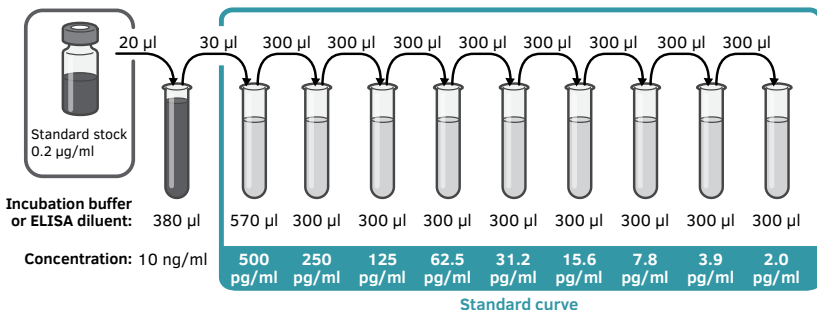


Protocol

Day 1

1. Add 100 μl /well of capture mAb 1-D1K diluted to 2 $\mu\text{g}/\text{ml}$ in PBS, pH 7.4. Use high protein binding ELISA plates. Incubate overnight at 4-8 $^{\circ}\text{C}$.

Day 2

2. Empty the plate and add 200 μl /well of PBS with 0.05% Tween 20 and 0.1% BSA (incubation buffer) to block the plate. Incubate for 1 hour at room temperature.
3. Wash the plate 5 times with PBS containing 0.05% Tween 20 (300 μl /well).
4. Add 100 μl /well of samples or standards diluted in incubation buffer or ELISA diluent. Include assay background control, i.e. wells without standard. Incubate for 2 hours at room temperature.
5. Wash as above.
6. Add 100 μl /well of detection mAb 7-B6-1-biotin diluted to 1 $\mu\text{g}/\text{ml}$ in incubation buffer or ELISA diluent. Incubate for 1 hour at room temperature.
7. Wash as above.
8. Add 100 μl /well of Streptavidin-HRP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature. Please note that sodium azide used in buffers will inhibit HRP activity.
9. Wash as above.
10. Add 100 μl /well of TMB substrate (product code: 3652-F10) and incubate at room temperature, protected from direct light for 15 minutes.
11. Add 100 μl /well of 0.2 M H_2SO_4 to stop the reaction.
12. Measure the optical density in an ELISA reader at 450 nm within 15 min. Preferably use a reader capable of subtracting a reference wavelength of between 570 and 650 nm. Representative standard curve shown below.

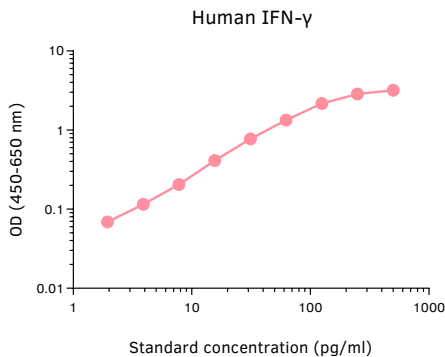

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.



The products are for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office)
Sweden
Tel: +46 8 716 27 00
mabtech@mabtech.com

Mabtech, Inc.
USA
Tel: +1 513 871-4500
mabtech.usa@mabtech.com