

ELISA Flex: Human IL-29 (ALP)

3570-1A-6 | 3570-1A-20

ELISA Flex kit for quantitative determination of native and recombinant human IL-29 (IFN- λ 1) in solution, e.g. cell supernatant.

The kit includes		3570-1A-6 for 6 plates	3570-1A-20 for 20 plates
Capture mAb:	MT2H7 (0.5 mg/ml)	300 μ l	1000 μ l
Detection mAb:	MT6G4, biotin (0.5 mg/ml)	150 μ l	500 μ l
Streptavidin-ALP		80 μ l	250 μ l
Recombinant human IL-29 (IFN- λ 1) ELISA standard		1 vial	1 vial
Standard reconstitution buffer A8		1 ml	1 ml

To ensure total recovery of the stated quantity, vials have been overfilled.

Shipping and storage

Shipped at ambient temperature. All reagents should be stored at 4-8 °C upon receipt, except the standard which should be stored at -20 °C. Antibodies are supplied in sterile-filtered PBS with sodium azide (0.02%). Streptavidin-ALP is supplied in 0.1 M Tris buffer with 0.002% Kathon CG. The expiry date indicates how long unopened products, stored according to instructions, are recommended for use.


Protocol

Day 1

1. Add 100 μ l/well of capture mAb MT2H7 diluted to 2 μ g/ml in PBS, pH 7.4. Use high protein binding ELISA plates. Incubate overnight at 4-8 $^{\circ}$ C.

Day 2

2. Empty the plate and add 200 μ l/well of PBS with 0.05% Tween 20 and 0.1% BSA (incubation buffer) to block the plate. Incubate for 1 hour at room temperature.
3. Wash the plate 5 times with PBS containing 0.05% Tween 20 (300 μ l/well).
4. Add 100 μ l/well of samples or standards diluted in incubation buffer or ELISA diluent. Include assay background control, i.e. wells without standard. Incubate for 2 hours at room temperature.
5. Wash as above.
6. Add 100 μ l/well of detection mAb MT6G4-biotin diluted to 1 μ g/ml in incubation buffer or ELISA diluent. Incubate for 1 hour at room temperature.
7. Wash as above.
8. Add 100 μ l/well of Streptavidin-ALP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature.
9. Wash as above.
10. Add 100 μ l/well of pNPP substrate (product code: 3652-P10) and incubate at room temperature protected from direct light for approximately 60 minutes.
11. Measure the optical density in an ELISA reader at 405 nm. Preferably use a reader capable of subtracting a reference wavelength of between 570 and 650 nm. Representative standard curve shown below.

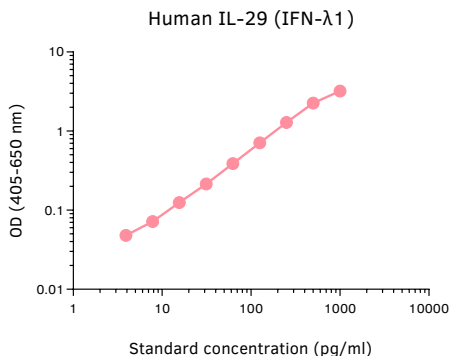

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.



The products are for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office)
Sweden
Tel: +46 8 716 27 00
mabtech@mabtech.com

Mabtech, Inc.
USA
Tel: +1 513 871-4500
mabtech.usa@mabtech.com